

INSTRUCTION  
MANUAL

**la parfait**  
COSMETICS

## TABLE OF CONTENTS

---

<b>WHAT'S INSIDE</b>	<b>1</b>
<b>DETAILED EXPLANATION</b>	<b>2</b>
<b>INSTRUCTIONS FOR USE AND STORAGE</b>	<b>3</b>
<b>COLOR BENEFITS</b>	<b>4-5</b>
<b>WARNINGS AND INFORMATION</b>	<b>6</b>

---

800 689 4243

730 w. Hallandale Beach Blvd.

Suite #105, Hallandale Beach Fl,

33431

**la parfait**  
COSMETICS

## WHAT'S INSIDE THE BOX

---



LED mask



Remote control



Instructional manual



USB connection cable




Charger



Support and retention strap

## MASK AND REMOTE INPUTS AND BUTTONS

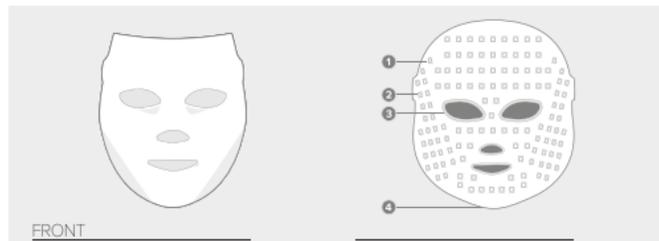

<b>1</b>	Mounted LED	Emits specific wavelengths and light colors.
<b>2</b>	Strap connectors	Used to secure the straps onto the mask.
<b>3</b>	Eye comforter	Rubber comfort pads for the eyes.
<b>4</b>	Power supply port	Power connection from the mask to the remote.



<b>A</b>	Charger cable port	Port that connects the remote to the power supply.
<b>B</b>	Power supply port	USB power connector to the mask.
<b>C</b>	Power LED light	LED light indicating when power is connected.
<b>D</b>	Power button	Button for turning the remote on and off.
<b>E</b>	Energy plus button	Increases the intensity of the LED lights.
<b>F</b>	Time minus button	Reduces the treatment time by ( 5 ) minutes.
<b>G</b>	Energy minus button	Reduces the strength of the LED lights
<b>H</b>	Screen	Shows treatment countdown and intensity.
<b>I</b>	Time LED light	Indicates that the timer has been activated.
<b>J</b>	ON/OFF button	Starts and stops the treatment.
<b>K</b>	Time plus button	Increases the treatment time by ( 5 ) minutes.
<b>L</b>	Color changing button	Changes the color/wavelength of the LED lights.

# INSTRUCTIONS FOR USE AND STORAGE

---

## MASK & TREATMENT SET UP

- 1 Connect the charger (see page 1) into **A** the charger cable port and plug the charger into the electrical wall socket. Press **D** the power button and check light **C** to see if it turned on. If the light is illuminated, you have successfully turned the device on.
- 2 Connect the USB connection cable (see page 1) into port **B** at the front of the remote and plug the opposite end in port **4** located at the chin of the led mask. The treatment is now ready to begin!
- 3 Press the ON/OFF button **J** to start the treatment. The timed LED light **I** will turn on and the screen **H** will show the countdown of the treatment in minutes.  
Find the correct LED color you desire by using **L** the button marked color. Once you found the correct treatment (see page 4) for information) you can adjust the LED light intensity level by using **G**, **E** and you can set the treatment time with **K** and **F**  
Depending on the skin condition. It is recommended to use everyday for five minutes or three times a week for 10 - 20 minutes.
- 4 When done with treatment press the ON/OFF button **J** and then the power button **D**. Unplug the mask, carefully disinfect it with wipes and let it dry before putting it away.

## CARE & STORAGE

- Store away in a dry, cool place.
- Clean and disinfect after every use.
- Handle with care.

## LED COLOR BENEFITS

### GREEN

Discoloration/neutralizing treatment - between red and blue light, is suitable for mixed skin. It improves circulation, which helps with cell healing, absorption and balance. It aids in evening out skin tone and improving discoloration and pigmentation issues resulting from environmental factors, such as sun damage. This treatment can help fade out the appearance of freckles, as well as of age spots, sun spots, reddish and brown patches, and liver spots.

### RED

Wrinkle reduction treatment - penetrates 8-10 mm deep into the skin and is suitable for any skin type, and especially beneficial for dull, dry, allergic skin. This treatment stimulates fiber cells, while the heat improves blood circulation and anti-oxidation and increases collagen production to reduce skin's pores and wrinkles and to make skin more elastic, smoother, and younger-looking. The mild frequency is suitable for the areas surrounding the eyes and forehead.

### BLUE

Spot/acne treatment - light penetrating about 0.5 mm has the function of calming inflamed skin to treat and prevent rashes. It is suitable for skin that is sensitive, oily and prone to rashes. It adjusts the skin's sebum secretion, effectively eliminating acne-causing bacteria and reducing acne and wheek in order to reach the ultimate goal of beautification. Additionally, the blue light astringes and tightens loosened skin.

### WHITE

Combination of all lights - should be used twice a week for maintenance of results. This treatment enhances skin's nutrient absorption, while helping minimize the appearance of wounds, acne, fine lines and wrinkles, and sagging skin.

## LED COLOR BENEFITS

---

### PURPLE

Combination of both blue light therapy for acne and red light therapy for wrinkles - the best option for maintenance and oilier, acne-prone skin, as well as for those concerned with persistent breakouts and the signs of aging. It increases skin's stem renewal and regeneration. Thus, this treatment treats acne and past acne marks. It also helps with lymphatic drainage..

### YELLOW

Energizes the skin cells and helps to soothe sensitive skin and brighten the complexion. This treatment also stimulates the lymphatic system, assists digestion, alleviates skin diseases, and enhances immune function.

### CYAN

Enhances cell energy intake and promotes metabolic functions. This treatment reduces skin tension, helps with healing of inflammation and sunburns, and relieves overall stress.

## WARNINGS AND INFORMATION

---

### OPERATING WARNINGS

- Use only with the power cables provided.
- Do not place any wet objects near the device.
- Do not operate with wet hands.
- Do not try to open or tamper with any part of the device.

### TREATMENT WARNINGS

- Do not use if pregnant
- Do not use whilst being treated for cancer
- Prior to use, seek medical advice in case of epilepsy or hypersensitivity to light.
- Do not look directly at the LED lights when in use.

<b>Product</b>	La parfait Led Face Mask
<b>Purpose / Intended usage</b>	To treat and help minimize symptoms of various skin conditions and problems.
<b>Quantity of Leds</b>	150 of each wavelength. Total of 1050 bulbs.
<b>Voltage and Power</b>	100 - 240VAC

LAPARFAITCOSMETICS.COM